
AIREX®

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AIREX®
CADEMY

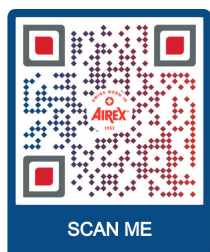AIREX®
YOGA & PILATES

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Sisällysluettelo / Table of Contents / Inhaltsverzeichnis

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	Valmistaja / Manufacturer Indicates the medical device manufacturer		Lämpötilan raja-arvot / Temperature limit Indicates the temperature limits to which the medical device can be safely exposed
	Valtuutettu edustaja / Authorized Representative in the EC Indicates the authorized representative in the European Community / European Union		Älä käytä, jos pakkaus on vaurioitunut / Do not use if package is damaged Indicates a medical device that shouldn't be used if the package has been damaged
	Valmistuspäivä / Date of manufacture Indicates the date when the medical device was manufactured		Noudata käyttöohjeita / Consult instructions for use Indicates the need for the user to consult the instructions for use
	Eränumero / Batch code Indicates the manufacturer's batch code so that the batch or lot can be identified		Varoitus / Caution To indicate that caution is necessary when operating the device or control close to where the symbol is placed
	Luettelonumero / Catalog number Indicates the manufacturer's catalog number so that the medical device can be identified		Lääkinnällinen laite / Medical device Indicates the item is a medical device
	Sarjanumero / Serial number Indicates the manufacturer's serial number so that a specific medical device can be identified		Käännös / Translation To identify that the original medical device information has undergone a translation
	Säilytä auringonvalolta suojattuna / Keep away from sunlight Indicates a medical device that needs protection from light sources		Tämä suunta ylöspäin / This way up
	Säilytä kuivassa / Keep dry Indicates a medical device that needs protection from moisture		CE-merkintä / CE marking 'CE marking of conformity' or 'CE marking' means a marking by which a manufacturer indicates that a device is in conformity with the applicable requirements
	Yksilöllinen laitetunniste / Unique device identifier Indicates a carrier that contains unique device identifier information		Maahantuoja / Importer To indicate the entity importing the medical device into the locale
			Jakelija / Distributor To indicate the entity distributing the medical device into the locale

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Great importance is attached to functionality in the super-soft products. A special manufacturing process is necessary to achieve good cushioning. Some deviations from perfection may occur with these products, such as occasional small surface bubbles, small deviations in the dimensional tolerance or slight discolorations. However, this does not affect the function in any way.

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To enjoy your gymnastics mat for a long time, we strongly recommend cleaning the mat at regular intervals, especially if it is used in extremely warm and humid environments (e.g., public, or private baths, thermal baths, wellness areas). This will prevent any changes in the surface structure and increase the service life. For cleaning, use a soft sponge or cloth with lukewarm soapy water. Please do not use scouring agents or scouring sponges. Avoid cleaning with petrol or other organic solvents. After cleaning, allow the mat to air dry thoroughly before storing.

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The closed-cell AIREX® soft foams do not absorb water and have a permanent load-bearing capacity in water. Attention: Do not use as a swimming or rescue aid! The Sanitized® finish effectively protects the gymnastics mat against unpleasant odours and microbial infestation. It preserves the material, improves hygiene and increases comfort.*

If disinfection is nevertheless necessary for special reasons, please use a commercially available water-based disinfectant. Disinfection or sterilisation with superheated steam is not permitted.

* The antimicrobial Sanitized® hygiene function integrated in our mats protects the material against odours, bacteria as well as mould. Contains a biocidal product with the active ingredient Folpet. Concentration: 0.1 - 0.2 %.

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Single mats can be rolled up without further ado. Make sure that the grooves are facing outwards. Storage on our special holders is ideal for storing several mats, provided your mats are equipped with AIREX® eyelets. Please always hang the mats on both eyelets. If you lose eyelets, you can replace them without any problem. Replacement eyelets can be obtained from our service department. (see contact).

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Think of the environment when disposing of packaging material and product and please enquire at your local administration about the possibilities of environmentally friendly disposal via appropriate collection points.
AIREX[®] gymnastic mats are made of flexible foam based on polyvinyl chloride cellular materials (PVC) and are expanded with a nitrogen-releasing chemical blowing agent. The material does not contain phthalates.

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We provide the prescribed legal warranty for AIREX[®] gymnastic mats, unless otherwise specified below.

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Claims for damage to the AIREX[®] gymnastics mat due to improper use are excluded.
Claims for damages - irrespective of the legal grounds, in particular due to breach of obligations arising from the contractual obligation and from unauthorised actions - are excluded. This shall not apply insofar as liability is mandatory by law, in particular:
• in cases of intent or gross negligence;
• in cases of injury to life, body or health;
• under the Product Liability Act or

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• in the event of any other breach
of material contractual obligations;
In this case, the claim for compen-
sation is limited to the typical dam-
age.

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Airex AG
3A Composites GmbH